
seca
Precision for health

seca 874

Instructions for Use

17-10-06-345-100g_2026-03S

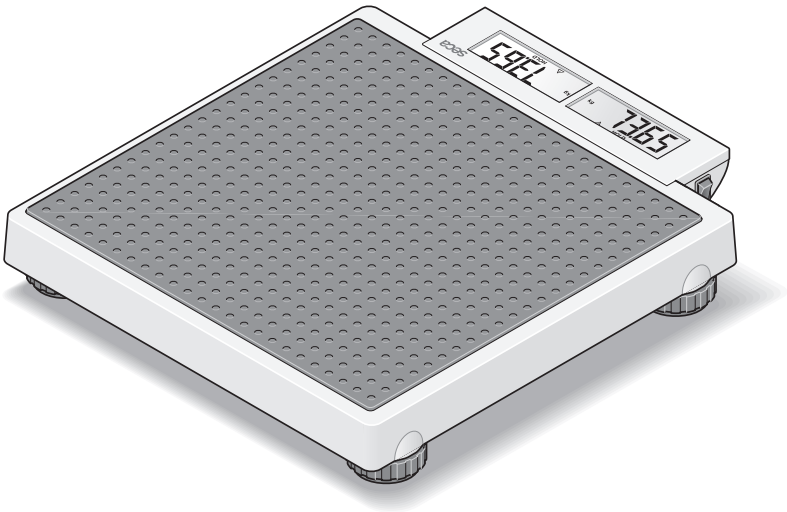seca.com

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English – Instructions for Use **27**

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17-10-06-345-100g_2026-03S

CONTENTS

1 Description of device 28
1.1 Intended use..... 28
1.2 User qualification 28
1.3 Patient target group..... 28
1.4 Contraindications..... 28
1.5 Clinical benefit 28
1.6 Functional description..... 29

2 Safety precautions..... 29
2.1 Safety information in these instructions for use 29
2.2 Basic safety information..... 30
Handling device 30
Preventing electric shock 31
Preventing injuries and infections 31
Preventing damage to device..... 32
Handling measuring results..... 33
Handling packaging material..... 34
Handling batteries and rechargeable batteries ... 34

3 Overview..... 35
3.1 Controls..... 35
3.2 Symbols in the display 36
3.3 Labels..... 36

4 Starting up device..... 39
4.1 Scope of delivery 39
4.2 Establishing power supply . 39
Inserting batteries 39
4.3 Setting up and aligning the device..... 40

5 Operation..... 40
5.1 Switching the device on and off 41
Switching on the device . 41
Switching the display on and off (standby)..... 41
Switching off the device . 41
5.2 Weighing..... 42
Weighing a patient 42
Weighing babies/ toddlers (2 in 1)..... 42
Switching the weight unit 43

6 Hygiene treatment 44
6.1 Cleaning 45
6.2 Disinfecting 45
6.3 Sterilizing 45

7 Function check..... 46

8 Troubleshooting..... 46

9 Servicing..... 48

10 Technical data..... 48
10.1 General technical data 48
10.2 Weight measurement..... 49

11 Disposal 49
11.1 Disposing of the device..... 49
11.2 Disposing of batteries and rechargeable batteries 49

12 Warranty 50

13 Declaration of conformity 50

17-10-06-345-100g_2026-03S

1 DESCRIPTION OF DEVICE

1.1 Intended use

The electronic flat scale supports physicians in decision-making regarding diagnosis or therapy, based on the patient's weight or weight changes.

1.2 User qualification

The device may only be used by persons with sufficient expertise.

1.3 Patient target group

The device is intended for persons of any age who can stand independently and who do not exceed the maximum weight capacity of the device. Babies and toddlers who cannot stand independently can be weighed using the **2 in 1** function whilst being held by an adult who can stand independently.

1.4 Contraindications

No contraindications are known.

1.5 Clinical benefit

The electronic flat scale supports physicians in decision-making regarding diagnosis or therapy, based on the measured weight.

By that, the electronic flat scale only provides an indirect clinical benefit.

The physician must perform additional examinations for an accurate diagnosis.

17-10-06-345-100g_2026-03S

1.6 Functional description

Weight on this electronic personal scale is detected by four load cells.

The **2 in 1** function can be used to determine the weight of babies and toddlers. The child is held in an adult's arms during weighing to do this.

2 SAFETY PRECAUTIONS

2.1 Safety information in these instructions for use

DANGER!
Used to identify an extremely hazardous situation. If you fail to take note of this information, serious irreversible or fatal injuries will occur.

WARNING!
Used to identify an extremely hazardous situation. If you fail to take note of this information, serious irreversible or fatal injuries may result.

CAUTION!
Used to identify a hazardous situation. If you fail to take note of this information, minor to moderate injuries may result.

NOTICE!
Used to identify possible incorrect usage of the device. If you fail to take note of this information, you may damage the device, or the measuring results may be incorrect.

NOTE
Includes additional information about use of the device.

17-10-06-345-100g_2026-03S

2.2 Basic safety information

Handling device

- ▶ Please take note of the information in these instructions for use.
- ▶ Keep the instructions for use in a safe place. The instructions for use are a component of the device and must be available at all times.
- ▶ In the interest of patient safety, you and your patients are obliged to report serious events that occur in connection with this product to the manufacturer and to the authority responsible in your country.

DANGER!
Risk of explosion

- ▶ Do not use the device in an environment in which one of the following gases has accumulated:
 - oxygen
 - flammable anesthetics
 - other flammable substances/air mixtures

CAUTION!
Patient hazard, damage to device

- ▶ Have servicing carried out regularly as described in the relevant section of this document.
- ▶ Technical modifications may not be made to the device. The device does not contain any user-serviceable parts. Only have servicing and repairs performed by an authorized seca service partner. You can find your local service partner at www.seca.com.
- ▶ Use only seca original accessories and spare parts, otherwise seca will not grant any warranty.

CAUTION!
Patient hazard, malfunction

- ▶ Keep medical electrical devices such as high-frequency surgical equipment at a minimum distance of approx. 1 meter to prevent faulty measurements.
- ▶ Keep high-frequency equipment such as cellphones at a minimum distance of approx. 1 meter to prevent faulty measurements.

17-10-06-345-100g_2026-03S

- The actual transmission output of HF devices may require minimum distances of more than 1 meter. For details, go to www.seca.com.

Preventing electric shock

WARNING!
Electric shock

- Set up devices which can be operated with a power supply unit so that the power supply socket is within easy reach and the power supply can be disconnected quickly.
- Ensure that your local electricity supply matches the details on the power supply unit.
- Never touch the power supply unit with wet hands.
- Do not use extension cables and multiple outlets.
- Ensure that cables are not crushed or damaged by sharp edges.
- Ensure that cables do not come into contact with hot objects.
- Do not operate the device at an altitude of more than 3000 m above sea level.

Preventing injuries and infections

WARNING!
Injury from falling

- Ensure that the device is steady and level.
- Route connector cables (if present) so that neither users nor the patient can trip over them.
- Assist people with limited motor skills when they are getting up, e.g. from a wheelchair.
- Ensure that the patient does not step directly onto or off the edges of the weighing platform.
- Ensure that the patient steps onto and off the weighing platform slowly and safely.

17-10-06-345-100g_2026-03S

WARNING!
Danger of slipping

- ► Ensure that the weighing platform is dry before the patient steps onto it.
- ► Ensure that the patient has dry feet before stepping onto the weighing platform.
- ► Ensure that the patient steps onto and off the weighing platform slowly and safely.

WARNING!
Risk of infection

- ► Before and after every measurement, wash your hands to reduce the risk of cross-contamination and nosocomial infections.
- ► Should the patient have contagious diseases, subject the device to a hygiene treatment as described in the relevant section of these instructions for use immediately afterwards.
- ► Ensure that the patient does not have any open wounds or infectious skin alterations which may come into contact with the device.
- ► Subject the device to a hygiene treatment at regular intervals.

Preventing damage to device

NOTICE!
Damage to device

- ► Ensure that fluids never get inside the device. These can destroy the electronics.
- ► For devices with power supply operation: Switch off the device before you disconnect the power supply unit from the power supply socket.
- ► For devices with power supply operation: If the device is not to be used for an extended period, disconnect the power supply unit from the power supply socket. Only then is the device de-energized.
- ► For devices with battery or rechargeable battery operation: If you are not using the device for an extended period of time, remove batteries or rechargeable batteries. Only then is the device de-energized.
- ► Do not drop the device.

17-10-06-345-100g_2026-03S

- ▶ Do not subject the device to shocks or vibrations.
- ▶ Perform a function check before each use as described in the corresponding section in this document. Do not operate the device if it is not working properly or is damaged.
- ▶ Do not place the device in direct sunlight and ensure that it is not placed in the direct proximity of a heat source. The excessive temperatures could damage the electronics.
- ▶ Avoid rapid temperature fluctuations. If the device is transported so that a temperature difference of over 20 °C occurs, the device must be left to stand for at least 2 hours before it is switched on, otherwise condensation may form; this may damage the electronics.
- ▶ Use the device only in the intended ambient conditions outlined in the section entitled “Technical data”.
- ▶ Store the device only in the intended storage conditions outlined in the section entitled “Technical data”.
- ▶ Use only cleaning agents and disinfectants which match the details in the section entitled “Hygiene treatment”.
- ▶ For scales: Ensure that maximum capacity is not exceeded.

Handling measuring results

NOTICE!
Inconsistent measuring results

- ▶ Before you save and continue using measured values determined using this device (e.g. in seca software or in an EMR system), ensure that the measured values are plausible.
- ▶ If measured values have been sent to seca software or to an EMR system, ensure before continuing to use them that the measured values are plausible and assigned to the correct patient.

17-10-06-345-100g_2026-03S

Handling packaging material

WARNING!
Risk of suffocation
Packaging materials made of plastic film (e.g. bags) present a risk of suffocation.

- ► Store packaging material out of the reach of children.
- ► If the original packaging material is no longer available, only use plastic bags with safety holes to reduce the risk of suffocation. Use recyclable materials if possible.

NOTE
Store the original packaging material for future use (e.g. returning for servicing).

Handling batteries and rechargeable batteries

WARNING!
Injury
Batteries and rechargeable batteries contain harmful substances which may explode and so escape if not handled properly.

- ► Do not try to recharge batteries.
- ► Do not expose batteries/rechargeable batteries to heat.
- ► Do not burn batteries/rechargeable batteries.
- ► If harmful substances are leaking out, avoid contact with the skin, eyes, and mucous membranes. Rinse affected areas of the body with plenty of clean water and seek medical help at once.

NOTICE!
Damage to device and malfunction due to incorrect handling


- ► Use only the type of battery/rechargeable battery quoted in this document.
- ► Always replace all batteries/rechargeable batteries at the same time.
- ► Do not short-circuit batteries/rechargeable batteries.
- ► If you do not use the device for an extended period of time, remove the batteries/rechargeable batteries. This prevents acid from leaking into the device.

17-10-06-345-100g_2026-03S

► If harmful substances have penetrated the device, do not continue using it. Have the device checked by an authorized seca service partner and repaired if necessary.

3 OVERVIEW


3.1 Controls

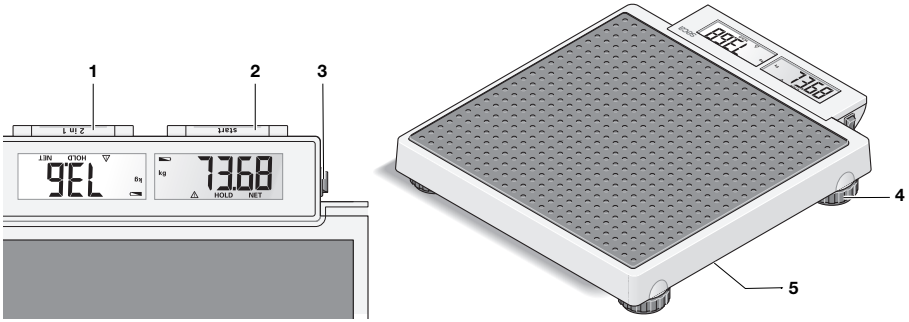Item	Control	Function
1	**2 in 1** key	Start the **2 in 1** function for weighing babies and toddlers
2	**start** key	Starts the weighing function
3	On/off switch	Switches the scale on and off
4	Foot screw	4 pcs, for precise alignment
5	Battery compartment	Takes 6 AA-type mignon batteries (AA/1.5 V)

17-10-06-345-100g_2026-03S

3.2 Symbols in the display

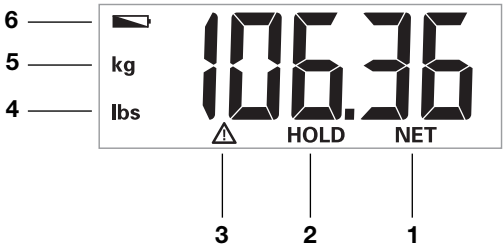Item	Symbol	Meaning
1	**NET**	**Tare** function active
2	**HOLD**	**Hold** function active
3	⚠	Non-verifiable function active
4	**lbs**	Weight value in pounds
5	**kg**	Weight value in kilograms
6		Batteries are low

3.3 Labels

Markings on the device and on the type plate	
Text/symbol	**Meaning**
	Name and address of manufacturer, date of manufacture
UDI	Unique Device Identification (product identification number as per Directive (EU) 2017/745)
REF	Article number
SN	Serial number
ProdID	seca product identification number
	Follow instructions for use

17-10-06-345-100g_2026-03S

Markings on the device and on the type plate	
Text/symbol	**Meaning**
	Medical electrical device, Type B
IP20	Type of protection to IEC 60529: • Protection against ingress of solid foreign bodies with a diameter of over 12.5 mm • Protection against access with fingers • No protection from ingress of water
d	Value in units of mass States the difference between two consecutive display values
CE 0123	Device complies with EU directives **0123**: Notified Body for Medical Devices
MD	Medical device in accordance with Regulation (EU) 2017/745
UK CA XXXX	Device complies with United Kingdom directives **xxxx**: Notified Body for Medical Devices of the United Kingdom
 UK REP	Importer/representative in the United Kingdom: seca Ltd 40 Barn Street B5 5QB Birmingham United Kingdom
 CH REP	Importer/representative in Switzerland: seca ag (schweiz) Medizinische Waagen und Messsysteme Schönmatt Str. 2 CH-4153 REINACH
	Do not dispose of device in household waste


17-10-06-345-100g_2026-03S

Markings on the packaging	
Symbol	**Meaning**
	Protect from moisture
	Arrows indicate top of product Transport and store in an upright position
	Fragile, do not throw or drop
	Permitted min. and max. temperature for transport and storage
	Permitted min. and max. humidity for transport and storage
	Permitted min. and max. air pressure for transport and storage
	Open packaging here
	Packaging material can be disposed of through recycling programs

17-10-06-345-100g_2026-03S

4 STARTING UP DEVICE


4.1 Scope of delivery

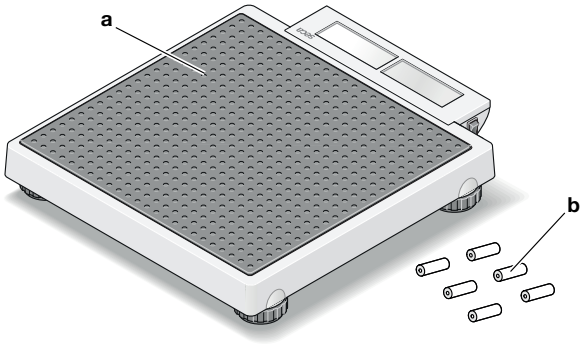Item	Component	Pcs.
a	Personal scale	1
b	Battery	6
–	Instructions for use	1

4.2 Establishing power supply

Inserting batteries

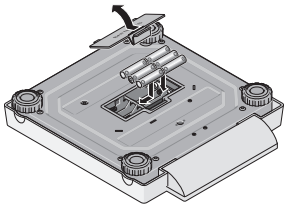You need 6 batteries (type AA/1.5 V, included in scope of delivery).

1. Press on the latch of the battery compartment.
2. Remove the cover of the battery compartment.

NOTICE!
Damage to device
Device may overheat and become damaged.

► Ensure the correct polarity of the batteries (markings in the battery holder).

3. Insert the batteries.

NOTE
If **bAtt** appears in the display, you have inserted one of the batteries the wrong way round or the batteries have discharged.


4. Close the battery compartment.

17-10-06-345-100g_2026-03S

4.3 Setting up and aligning the device

CAUTION!
Faulty measurement as a result of devices being set up incorrectly
Soft floors (wooden boards, for example) give under the patient's weight and falsify the measuring result. Incorrect ambient conditions or lack of alignment lead to measurement errors.

- ► Select a setup location with a level, stable floor in order to achieve precise measuring results.
- ► Use the device only in the intended ambient conditions outlined in the section entitled “Technical data”.
- ► Align the device before using it and every time it changes location.

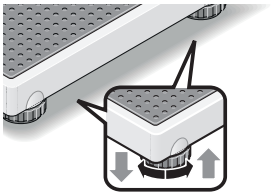1. Place the device on a firm, level surface.
2. Align the device by turning the foot screws.
3. Ensure that only the feet of the device are in contact with the ground.

5 OPERATION

CAUTION!
Personal injury/faulty measurement

- ► Before using the device each time, perform a function check as described in the corresponding section of these instructions for use.


17-10-06-345-100g_2026-03S

5.1 Switching the device on and off

Switching on the device

✓ The scale has no load.

► Set the ON/OFF switch to position “I”.

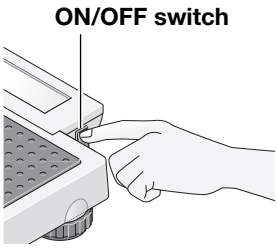
⇨ **SECA** appears in the display.

⇨ All the elements of the display are shown briefly.

⇨ The scale is operational when **0.00** appears in the display.

Switching the display on and off (standby)

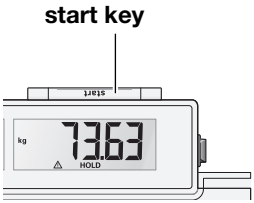
► Press the **start** key.

NOTE
The display switches off automatically after a short time if there is no load on the scale (not when operated with a power supply unit).

Switching off the device

► Set the ON/OFF switch to position “O”.

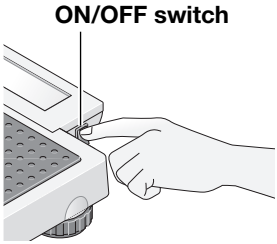NOTE
Switch off the device before transporting it so that the **start** key is not accidentally pressed during transport, so wasting power.

17-10-06-345-100g_2026-03S

5.2 Weighing


Weighing a patient

✓ The scale has no load.

1. If the scale switched off automatically after the last measurement, press the **start** key.

NOTICE!
Damage to device
The display is sensitive and may be damaged under load

► Make sure that the patient does not step on the display.

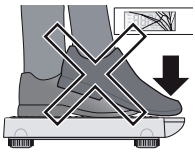2. Ask the patient to step onto the scale.

NOTE
The HOLD function is activated automatically at weights of over 1.5 kg or 3.3 lbs. The display flashes until a stable weight has been measured. The display then freezes until the next weighing operation. ⚠ **HOLD** appears in the display.


3. Read off the measuring result.

Weighing babies/ toddlers (2 in 1)

The **2 in 1** function can be used to determine the weight of babies and toddlers. The child is held in an adult's arms during weighing to do this. Proceed as follows:

1. Ensure that there is no load on the scale.
2. Ask the adult to stand on the scale.
3. Read off the measuring result.

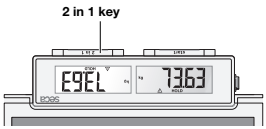4. Keep the **2 in 1** key depressed until the message **NET** appears in the display.
 ⇨ The weight is saved.
 ⇨ **0.00** appears in the display.

17-10-06-345-100g_2026-03S

WARNING!
Faulty measurement as a result of initial weight changing
If the child is weighed together with a different adult, it will not be possible to determine the child's weight correctly.

- ► Ensure that the child is always weighed with the adult with whom the initial weight was determined.
- ► Ensure that the adult's weight does not change - as a result of their taking off a piece of clothing, for example.

5. Ask the adult to get onto the scale holding the child.
 ⇨ The child's weight is displayed.
 ⇨ The symbol ⚠ and the messages **HOLD** and **NET** are displayed.
6. Read off the measuring result.
7. Ask the adult to get off the scale with the child.
8. You have the following options for deactivating the **2 in 1** function:
 - ► Press the **2 in 1** key until the symbol ⚠ and the messages **HOLD** and **NET** are no longer displayed.
 - ► Switch off the device.


Switching the weight unit

You can select the unit in which weight is to be displayed.

WARNING!
Patient hazard
To prevent misinterpretations, measuring results for medical purposes must only be displayed and used in SI units (kilograms/grams, meters/centimeters).


- ► Only use the measuring results in SI units.
- ► You as the user take sole responsibility for the use of measuring results in non-SI units.

17-10-06-345-100g_2026-03S

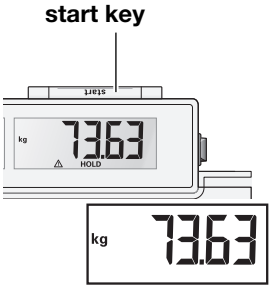1. Press the **start** key.

⇨ The current weight unit (in this case kg) will be displayed.

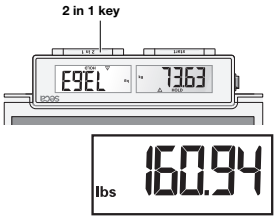2. Keep the **2 in 1** key depressed for approx. 3 seconds.

⇨ The unit is switched (in this case lbs).

3. To switch the unit back again, keep the **2 in 1** key depressed for approx. 3 seconds again.

6 HYGIENE TREATMENT

WARNING!
Electric shock
The device is not de-energized when the on/off button is pressed and the display goes out. Use of fluids on the device may cause an electric shock.

- ► Before each hygiene treatment, ensure that the device is switched off.
- ► Take the batteries/battery packs out of the device before every hygiene treatment.
- ► Ensure that no fluids penetrate the device.

WARNING!
Risk of infection

- ► Subject the device to a hygiene treatment at regular intervals as described in this section.

17-10-06-345-100g_2026-03S

NOTICE!
Damage to device
Unsuitable cleaning agents and disinfectants may damage the sensitive surfaces of the device and impair its operability.

- ► Do not use aggressive or abrasive cleaning agents.
- ► Do not use organic solvents (e.g. white spirit or petroleum spirit).

6.1 Cleaning

- ► If required, moisten a soft cloth with a mild soap solution and wipe the device over with it.

6.2 Disinfecting

1. Disinfect the device at regular intervals using a disinfectant suitable for sensitive surfaces and acrylic glass. seca recommends a 70 % ethanol solution as a surface disinfectant.
2. Follow the instructions for use of the disinfectant.
3. Disinfect the device as follows:
 - ► Moisten a soft cloth with disinfectant and wipe down the device with it.
 - ► Comply with the intervals, see table:

Interval	Component
Before and **after** every measurement with direct skin contact	Weighing platform
As required	Display

6.3 Sterilizing

The device must not be sterilized.

17-10-06-345-100g_2026-03S

7 FUNCTION CHECK

► Perform a function check before each use.

A complete function check includes:

- A visual inspection for mechanical damage
- A test of the alignment of the device
- A visual and function check of the display elements
- A function check of all the controls shown in the section entitled “Overview”
- A function check of the optional accessories

If you find faults or deviations during the function check, first try to remedy the fault with the aid of the “Troubleshooting” section in this document.

CAUTION!
Personal injury
If you find faults or deviations during the function check which you are unable to remedy with the aid of the “Troubleshooting” section in this document, you must not use the device.

► Have the device repaired by seca Service or an authorized service partner.

► Follow the section entitled “Servicing” in this document.

8 TROUBLESHOOTING

If faults occur when operating the device, first attempt to remedy them yourself using the following table. If the fault persists, contact seca Service or an authorized service partner.

17-10-06-345-100g_2026-03S

Fault	Cause	Remedy
The scale does not display a weight although a load was placed on it	Scale is switched off or without batteries	• Check that the scale is switched on • Check that batteries are inserted
Scale switches itself on, e.g. during transport	**start** key pressed	Set the On/off switch to “O” when transporting
The scale displays a weight after transport or a battery replacement although no load has been placed on it	Scale has determined an incorrect zero point	Switch off the scale, then switch it back on again
0.00 does not appear before weighing	Scale had load on it before being switched on	• Remove the load from the scale • Switch off the scale, then switch it back on again
---- displayed before weighing	Scale had load on it before being switched on	• Remove the load from the scale • Switch off the scale, then switch it back on again
Single segment is illuminated continuously or not at all	Display faulty	Inform seca Service
Display: 	Battery voltage is dropping	Replace the batteries
Display: **bAtt**	Batteries discharged or not inserted	• Replace the batteries • Insert the batteries the right way round
Display: **StOP**	Maximum capacity exceeded	Remove the load from the scale
Display: **Er** + figure	Too high a load on the scale or too high a load in one corner	• Remove the load from the scale or distribute the weight more evenly • Switch off the scale, then switch it back on again • If the problem persists, inform seca Service

17-10-06-345-100g_2026-03S

9 SERVICING

The product needs to be set up carefully and serviced regularly. Depending on how frequently the product is used, seca recommends servicing at intervals of 3 to 5 years.

CAUTION!
Faulty measurements as a result of poor servicing

- ▶ Have servicing and repairs carried out exclusively by seca Service or an authorized service partner.
- ▶ You can find your local service partner at www.seca.com.

10 TECHNICAL DATA

10.1 General technical data

Dimensions • Depth • Width • Height	 360 mm 321 mm 60 mm
Net weight	approx. 4.0 kg
Ambient conditions, operation • Temperature • Air pressure • Humidity	 +10 °C to +40 °C 700 – 1060 hPa 30 % – 85 %, no condensation
Ambient conditions, storage/transport • Temperature • Air pressure • Humidity	 -10 °C to +65 °C 700 – 1060 hPa 0 % – 95 %, no condensation
Height of digits	20 mm
Power supply • Batteries – Supply voltage – Battery type	 9 V 6 x type AA, 1.5 V, not rechargeable
Maximum current consumption	Typically 15 mA
Medical device in accordance with Regulation (EU) 2017/745	Class I with measuring function

17-10-06-345-100g_2026-03S

EN 60601-1: Medical electrical equipment, type B	
Application part in line with IEC 60601-1	Tread area
Type of protection	IP20
Duty cycle	Continuous duty

10.2 Weight measurement

Maximum capacity	200 kg/440 lbs
Accuracy • 0 to 50 kg/110 lbs • 50 kg/110 lbs to 200 kg/440 lbs	 ± 100 g/0.22 lbs ± 0.15 %
Graduation • 0 to 150 kg/330 lbs • 150 kg/330 lbs to 200 kg/440 lbs	 50 g/0.1 lbs 100 g/0.2 lbs

11 DISPOSAL

11.1 Disposing of the device

Do not dispose of the device in your household waste. The device must be properly disposed of as electronic scrap. Follow your respective national regulations. For more information, please contact seca Service or an authorized service partner.

11.2 Disposing of batteries and rechargeable batteries

Do not dispose of used batteries and rechargeable batteries in household waste, regardless of whether they contain harmful substances or not. As a consumer, you are legally obliged to dispose of batteries and rechargeable batteries via local authority collection points or trade collection points. Only dispose of batteries and rechargeable batteries once they are completely discharged.

17-10-06-345-100g_2026-03S

12 WARRANTY

There is a two-year warranty period from delivery for defects attributable to poor materials or workmanship. All movable parts, e.g. batteries, cables, power supply units, rechargeable batteries etc. are exempt. Defects which come under the warranty will be repaired for the customer free of charge against proof of purchase. Additional claims cannot be considered. Costs of transport to and from seca are the responsibility of the customer if the device is located somewhere other than the customer's headquarters. In the event of transport damage, claims under warranty can only be made if the complete original packaging was used for transport and the device was secured and fastened in it according to its originally packaged condition. You should therefore keep all packaging parts.

The warranty will be voided if the device is opened by persons not expressly authorized by seca to do so.

In the event of a warranty issue, please contact your local seca office or the dealer from whom you ordered the product.

13 DECLARATION OF CONFORMITY

CE 0123

seca gmbh & co. kg hereby declares that the product complies with the terms of the applicable European directives and regulations. The unabridged declaration of conformity can be found at www.seca.com.

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17-10-06-345-100g_2026-03S

17-10-06-345-100g_2026-03S

Medical Measuring Systems and Scales since 1840

seca gmbh & co. kg
Hammer Steindamm 3–25
22089 Hamburg · Germany
T +49 40 20 00 00 0
F +49 40 20 00 00 50
E info@seca.com

seca operates worldwide with headquarters
in Germany and branches in
seca france
seca united kingdom
seca north america
seca schweiz
seca zhong guo
seca nihon
seca mexico
seca austria
seca polska
seca middle east
seca suomi
seca américa latina
seca asia pacific
seca danmark
seca benelux
seca lietuva
and with exclusive partners in
more than 110 countries.
All contact data at seca.com